



FOLEY BALLOON CATHETER

MW 500 (2 WAY) | MW 501 (3 WAY)

- Foley Balloon Catheter is a Urinary tract catheter used for short or long term urine drainage.
- Catheter are made from Natural Latex Rubber.
- Siliconised smooth surface for non traumatic intubation.
- Distal end is coned with super smooth eyes proportional to each other to provide non traumatic intubation.
- Provided with ultra thin, highly elastic Balloon.
- Hard Non Return Valve ensures trouble free inflation and deflation of Balloon.
- Also available in 100% Silicon elastomer.
- For sizes 16 to 26 Fr, smaller balloon capacity is also available on customer requirement.



Packing

- Sterile, Individually packed in peelable blister pack.
- Box of 10; Master Carton of 600
- Gross weight: 15.300kg
- Master Carton size: 680 x 315 x 530mm

TWO WAY		THREE WAY	
Size in FG	Balloon Capty.	Size in FG	Balloon Capty.
6 & 10	3cc	16, 18, 20, 22, 24 & 26	30-50cc
12 & 14	15cc		
16, 18, 20, 22, 24 & 26	30-50cc		


SIZES										
FG	8	10	12	14	16	18	20	22	24	26
Colour Code	Black	Grey	White	Green	Orange	Red	Yellow	Voilet	Blue	Pink

NELATON CATHETER | MW 502

- Nelaton Catheter are used for bladder catheterisation through urethra.
- Catheter are made from Non toxic, non irritant, Medical grade PVC Compound.
- Radio opaque line provided through out the catheter for X-ray visualization.
- Distal end is coned with two lateral eyes for maximum drainage.
- Proximal end is provided with Funnel shape connector for easy connection to Urine Bag.
- Funnel shape connectors are colour coded for identification of sizes.

Packing

- Sterile, Individually packed in peelable blister pack.
- Box of 100; Master Carton of 1000.
- Gross weight: 14.200kg
- Master Carton size: 680 x 315 x 530mm

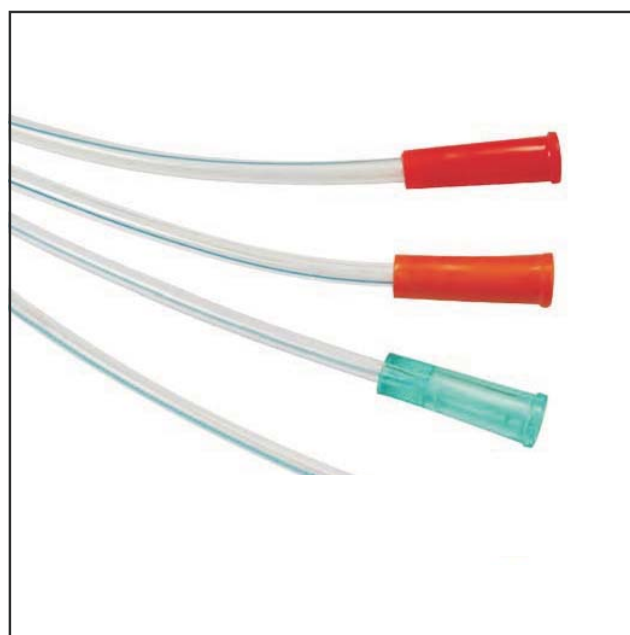

Versions	Effective Length	Total Length	Versions	Effective Length	Total Length
Male	370mm	400mm	Female	170mm	200mm

SIZES								
FG	6	8	10	12	14	16	18	20
Colour Code	Green	Blue	Black	White	Green	Orange	Red	Yellow



URINE COLLECTING BAG WITH MEASURED VOLUME METER | MW 510

- Urine Collecting Bag with measured volume meter is used for measurement of urine output.
- 2000ml Urine Bag is manufactured from extra strong sheeting to withstand during the long period.
- 250 ml Capacity measured Bottle provides accurate measurement of urine output.
- Urine Bottle is directly connected to the Urine Bag as to avoids any possibility of contamination.
- Push/Pull type Bottom outlet facilitates convenient emptying of Urine Bag.
- 90cm kink resistance tube provided with Universal connector.

Packing

- Sterile, Individually packed in peelable pouch pack.
- Box of 1; Master Carton of 50
- Gross weight: 5.500kg
- Master Carton size: 710 x 420 x 600mm



SIZES	
Bottle Capacity	250ml (Adult)
Bag Capacity	2000ml
Tube Length	90cm



URINE COLLECTING BAG | MW 504

- Urine Collecting Bag is used for short and long term urine drainage.
- 2000ml Urine Bag is manufactured from Medical grade soft PVC sheeting.
- Non Return Valve to avoid the back flow of Urine.
- Super smooth, 90cm kink resistance tube provided with Universal Connector for easy connection to the drainage catheter.
- Bag is graduated with easy to read scale.
- Also available with Top Bottom Outlet and Moulded handle for easy handling and carrying facility.

Packing

- Sterile, Individually packed in Tear H.M pack.
- Master Carton of 400.
- Gross weight: 25.200kg
- Master Carton size: 750 x 305 x 650mm

SIZES	
Bag Capacity	2000ml
Tube Length	90cm




URINE COLLECTING BAG PAEDIATRIC | MW 507

- Urine Collecting Bag Paediatric is used for urine collection for both Male and Female infants.
- 100ml Urine Bag consist of Hypo allergic adhesive for proper fixing of urine bag.
- Urine Bag is graduated with easy to read scale.

Packing

- Sterile, Individually packed in polybag.
- Box of 100; Master Carton of 2000.
- Gross weight: 18.100gm
- Master Carton size: 650 x 380 x 430mm

SIZES	
Bag Capacity	100ml

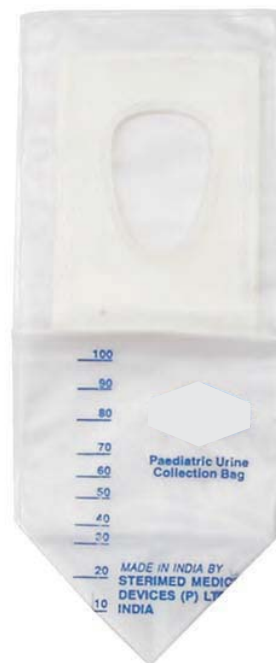

LEG BAG | MW 516

- Suitable for day & night use for incontinence
- Manufactured from soft clinical grade PVC sheeting.
- Available with bottom outlet.
- Non-return valve to prevent the backflow of urine.
- Printed with easy to read scale.
- Convenient for tying on thigh with optional elastic straps

Packing

- Sterile, Individually packed in peelable pouch pack.
- Box of 10; Master Carton of 100
- Gross weight: 9.200kg
- Master Carton size: 570 x 320 x 340mm

Capacity	Available Tube Lengths
750ml	10, 20 & 30cm
500ml	10, 20 & 30cm

LEVIN'S TUBE STOMACH TUBE | MW 602



- Levin's tube is used for Gastro intestinal drainage.
- Manufactured from non toxic, medical grade PVC compound.
- Distal end is coned with four lateral eyes for efficient drainage.
- The tube is marked at 45, 55, 65 and 75cm from the distal end for accurate placement into the abdomen.
- Radio opaque line provided throughout the tube for checking the position of the catheter radiologically.
- Proximal end is provided with Universal Funnel Connector for easy extension.

Packing

- Sterile, Individually packed in peelable pouch pack.
- Box of 50; Master Carton of 500
- Gross weight: 16.500kg
- Master Carton size: 790 x 475 x 275mm

SIZES								
FG	6	8	10	12	14	16	18	20
Colour Code	Green	Blue	Black	White	Green	Orange	Red	Yellow

Effective Length	119cm
Total Length	122cm

STERI-LENE

RYLE'S TUBE | MW 601

WITH LUER MOUNT | MW 601 L

- Ryle's tube is used for nutritional feeding purposes or aspiration of intestinal secretions.
- Manufactured from non toxic, medical grade PVC compound.
- Distal end is coned with steel balls sealed into tube to facilitates easy insertion.
- Four lateral eyes for efficient drainage.
- The tube is marked at 50, 60 and 70cm from the distal end for accurate placement into the abdomen.
- Radio opaque line provided throughout the tube for X-ray visualization.
- Proximal end is provided with Universal Funnel Connector for easy extension.
- Also available SMD 601 L with Luer mount fitting.

Packing

- Sterile, Individually packed in peelable pouch pack.
- Box of 50; Master Carton of 500
- Gross weight: 14.500kg
- Master Carton size: 790 x 275 x 475mm



SIZES								
FG	6	8	10	12	14	16	18	20
Colour Code	Green	Blue	Black	White	Green	Orange	Red	Yellow

Effective Length	102cm
Total Length	105cm



RECTAL CATHETER | MW 605

- Used for introduction of enema solution into rectum.
- Manufactured from non-toxic, non irritant medical grade PVC compound
- Distal end is open with one lateral eye for efficient drainage.
- Proximal end is provided with Funnel shape connector for easy connection to the suction equipment.
- Funnel connectors are colour coded for instant size identification.
- Also available in 40cm length.



Packing

- Sterile, Individually packed in peelable pack.
- Box of 100; Master Carton of 1000.
- Gross weight: 25.500kg
- Master Carton size: 750 x 305 x 650mm

SIZES								
FG	12	14	16	18	20	22	24	26
Colour Code	White	Green	Orange	Red	Yellow	Voilet	Blue	Pink

Effective Length	57cm
Total Length	60cm

ENDOTRACHEAL TUBE CUFFED

NASAL/ORAL | MW 701 C

- Manufactured from siliconised non-toxic, non irritant PVC kink resistant thermo sensitive tube ensures tube potency for patient safety. Yet softness at body temperature to conform to the anatomy of the respiratory tract.
- Murphy eye has a polished smoothness ensuring a traumatic intubation and extubation
- Kink resistant inflation tube ensures patient safety during cuff inflation and deflation.
- Radio opaque line and markings provided on the tube facilitate the exact location of the tube position
- Standard 15mm connector provided at proximal end E.T. Sterile & Pyrogen free

Packing

- Sterile, Individually packed in peelable pouch pack.
- Box of 10; Master Carton of 100.
- Gross weight: 4.500kg
- Master Carton size: 570 x 320 x 340mm



SIZES

2.5mm, 3.0mm, 3.5mm, 4.0mm, 4.5mm, 5.0mm, 5.5mm, 6.0mm, 6.5mm, 7.0mm, 7.5mm, 8.0mm, 8.5mm, 9.0mm, 9.5mm, 10mm

ENDOTRACHEAL TUBE PLAIN

NASAL/ORAL | MW 701 P

- Endotracheal Tube is used to provide a direct and unimpeded airway for gases to pass to and from the lungs.
- Manufactured from non toxic, non irritant medical grade PVC Compound.
- Tip is suitable for both Nasal and Oral intubation.
- Radio opaque line provided throughout the length for X-ray visualization.
- Proximal end is fitted with standard 15mm connector.

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 10; Master Carton of 100.
- Gross weight: 2.300kg
- Master Carton size: 570 x 320 x 340mm



SIZES


2.0mm, 2.5mm, 3.0mm, 3.5mm, 4.0mm, 4.5mm, 5.0mm, 5.5mm, 6.0mm, 6.5mm, 7.0mm, 7.5mm, 8.0mm, 8.5mm, 9.0mm, 9.5mm, 10mm

NASAL OXYGEN CANNULA TWIN BORE OXYGEN DELIVERY SET MW 703 A

- Nasal Oxygen Cannula is used for efficient administration of oxygen for greater comfort of patient.
- Twin Nasal Tips provide maximum freedom to the patient and leaves the mouth free.
- Nasal Tips ensures equal volume of oxygen to both air passage.
- Manufactured from soft non toxic, medical grade PVC Compound.
- Multichannel tube ensure regular supply of oxygen if tube is accidentally kinked.
- Proximal end of tube provided with connector for easy connection with oxygen source.

Packing

- Non-Sterile, individually packed in a polybag.
- Master Carton of 100.
- Gross weight: 6.00kg
- Master Carton size: 520 x 240 x 230mm

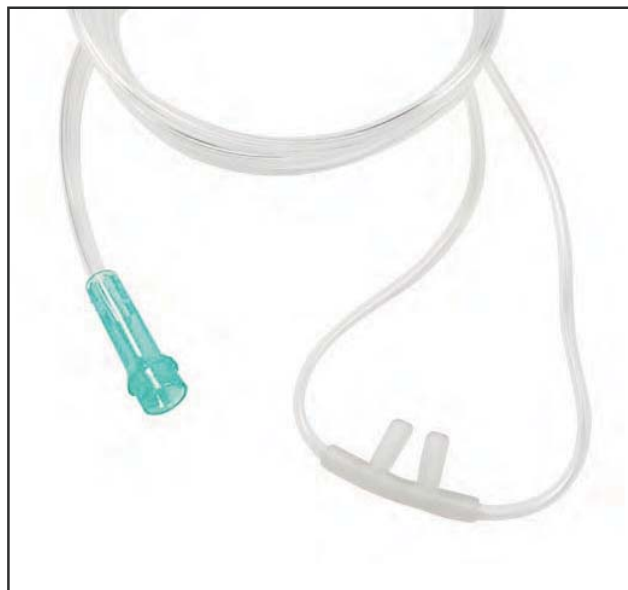

SIZES	
Size	Tube Length
Adult	210cm (7ft)
Child	210cm (7ft)
Infant	210cm (7ft)



OXYGEN MASK | MW 704

- Oxygen Mask is used for convenient oxygen administration.
- Manufactured from soft non toxic, medical grade PVC Compound.
- Elastic strap and Aluminum nose clip is provided on the Mask for proper adjustment of Mask on nasal area.
- Multichannel tube ensure regular supply of oxygen if tube is accidentally kinked.
- Tube is provided with connectors at both ends for easy connection with oxygen source.

Packing

- Non- Sterile, individually packed in a polybag.
- Master Carton of 50.
- Gross weight: 5.200kg
- Master Carton size: 520 x 240 x 230mm

SIZES	
Size	Tube Length
Adult	210cm (7ft)
Child	210cm (7ft)

SUCTION CATHETERS

MW 700 P (Plain)

MW 700A TC (Thumb Control)

MW 700B FT (Finger Tip Control)


MW 700 T (Automatic)

- Suction Catheters are used for removal of secretion from mouth, trachea and bronchial tube
- Manufactured from non toxic, non irritant medical grade PVC Compound.
- Distal end is open with one lateral eye for effective drainage
- Proximal end is provided with connector for easy connection to the suction equipment.
- Funnel connector are colour coded for instant size identification.
- Provided with vacuum control facility in Finger Tip Control and Thumb Control models


Packing

- Sterile, individually packed in peelable pack
- Box of 100; Master Carton of 1000.
- Gross weight: 18.500kg
- Master Carton size: 680 x 305 x 650mm

Versions	Efective Length	Total Length	Versions	Efective Length	Total Length
Plain	500mm	530mm	FT	470mm	530mm
TC	500mm	530mm	Automatic	500mm	530mm

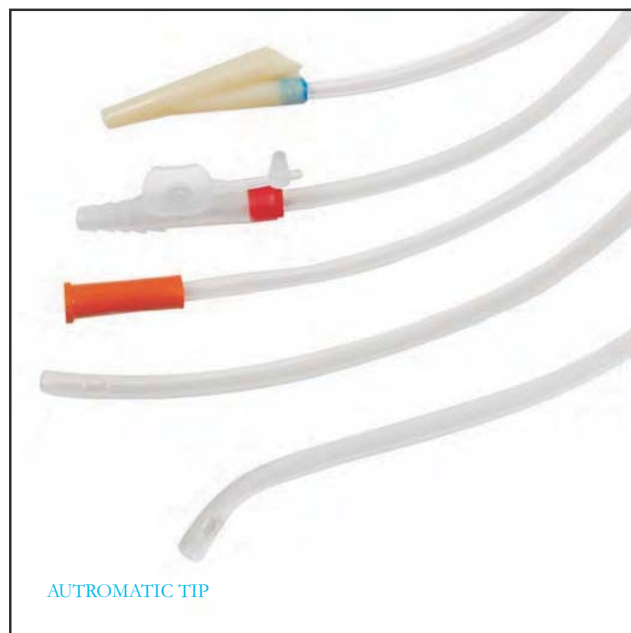

SIZES								
FG	5	6	8	10	12	14	16	18
Colour Code	Grey	Green	Blue	Black	White	Green	Orange	Red

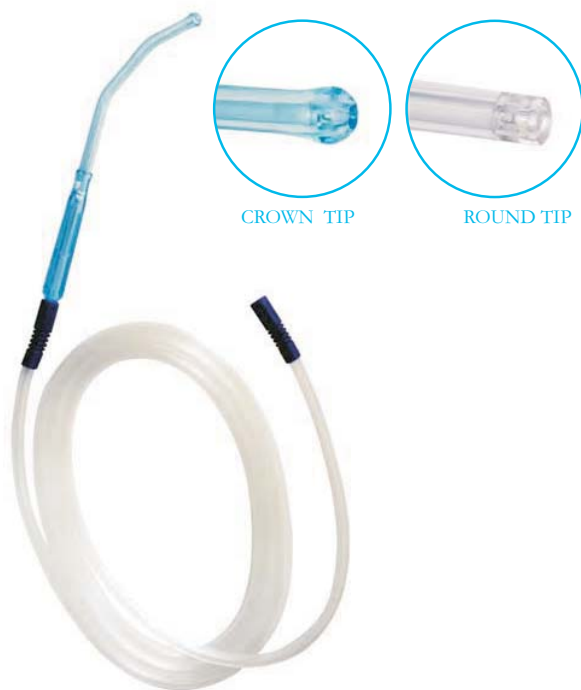

YANKAUR SUCTION SET

WITH STANDARD HANDLE MW 712

- Yankaur Suction Set is used for removal of secretion, and blood per operatively .
- Flexible kink resistant tube provide uninterrupted flow during suction.
- Soft connectors at both ends of the tube for easy connection.
- Available in Crown Tip Handle and Standard Round Tip Handle
- Suitable for use with or without vent control facility.

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 10; Master Carton of 100
- Gross weight: 19.200kg
- Master Carton size: 680 x 470 x 520mm

SIZES	
Tube Length	2500mm (2.5mtr.)
Tube Diameter	9.00mm

3 WAY STOP COCK | MW 708

- Three way stop cock are used for pressure infusion and blood pressure monitoring.
- Three arm handle for facilitating smooth and easy operations.
- Minimum residual fluid volume.
- Stop cock having two female luer lock ports and one male luer lock port with rotatory lock for safe connection to other device.
- 3 way stop cock available with extension line also.

Packing

- Box of 50; Master Carton of 500
- Gross weight: 6.200kg
- Master Carton size: 460 x 250 x 540mm



SIZES	
Chamber Capacity	25ml
Catheter Size	10 FG

INFANT MUCUS EXTRACTOR | MW 803

- Infant Mucus Extractor is used for Aspiration of secretion from oropharynx in newly born babies.
- Manufactured from non toxic, medical grade polymer.
- Transparent Chamber for visual examination of the content.
- Catheter is provided with round open tip for non traumatic intubation.
- Proximal end of suction tube is provided with Funnel shape connector for easy connection to the suction equipment.
- Also available with Filter and without Filter Models.

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 20; Master Carton of 400.
- Gross weight: 13.500kg
- Master Carton size: 770 x 410 x 490mm

GUEDEL AIRWAYS | MW 705

- Guedel Airways is used to maintain an unobstructed oro-pharyngeal airway.
- Bite blocks are fitted in the Airway to avoid airway occlusion and biting of tongue.
- Bite Blocks are colour coded for size identification.

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 50; Master Carton of 400.
- Gross weight: 9.500kg
- Master Carton size: 780 x 310 x 460mm




SIZES			
Size	Colour Code	Size	Colour Code
000	-	2	Green
00	Blue	3	Yellow
0	Black	4	Red
1	White	5	Orange

GYNAECOLOGY/PAEDIATRIC

UMBILICAL CATHETER | MW 1103

- Umbilical Catheter is used to access the umbilical vessels of pre-matured infants or newly born babies.
- Manufactured from non toxic, medical grade PVC.
- Smooth round tip for non traumatic insertion and well finished surface facilitates smooth passage of catheter in the vein.
- Radio opaque line provided throughout the tube for checking the position of the catheter radiologically.
- Marking is provided at every cm to ascertain depth of insertion.
- Proximal end of catheter is provided with Luer Female connector for easy connection to the I.V. Therapy equipment.

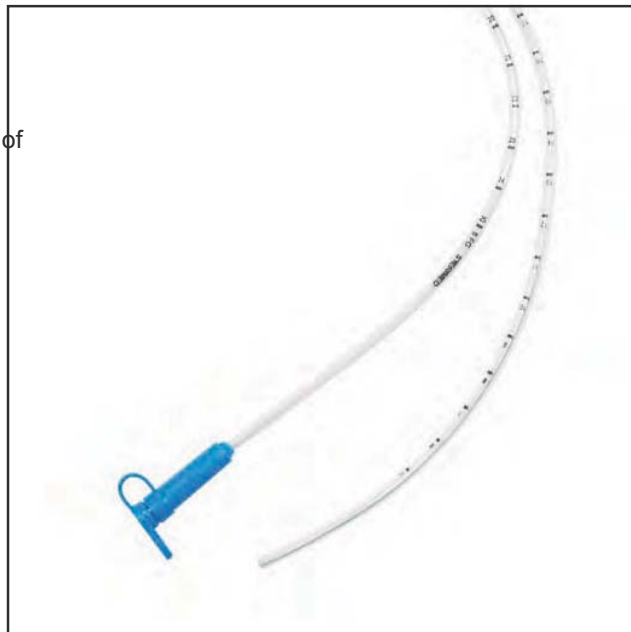

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 100; Master Carton of 1000
- Gross weight: 8.500kg
- Master Carton size: 680 x 310 x 530mm

Total Length - 35cm

SIZES					
FG	4	5	6	8	10
Colour Code	White	Grey	Green	Blue	Black



UMBILICAL CORD CLAMP | MW 801

- Umbilical Cord Clamp is used for clamping the umbilical cord of newly born baby.
- Manufactured from non toxic, medical grade polymer.
- Permanent locking system protect against the reopening of the clamp.
- Finger grip is provided to facilitate easy holding and locking of the clamp.
- Catheter clamp is also available for clamping the catheters.

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 100; Master Carton of 2000.
- Gross weight: 10.500kg
- Master Carton size: 650 x 380 x 430mm



STERIDERM I.V | MW116

- I.V Cannula Fixation Dressing
- Available in two variants: Transparent and Non Woven
- Two adhesive systems
- Reinforced notches reduce mechanical stress from heavy catheters
- Sterile tape strips enhance secural
- Polyurethane film is thin for better contour over high-profile devices
- Waterproof, sterile dressing cannot become saturated during bathing or showering
- Transparent, breathable and impermeable to liquids, viruses and bacteria

Packing

- Box of 50; Carton of 1000.
- Gross weight: 2.550kg
- Master Carton size: 220 x 150 x 220mm


THORACIC TROCAR CATHETER

MW 1001

- Thoracic Trocar Catheter is used for rapid non operative plueral and Chest drainage.
- Manufactured from non toxic, medical grade PVC compound.
- Catheter is snugly fxed over the trocar.
- Smooth taper shape Tip.
- Catheter is provided with large smooth eyes for eficient drainage.
- Catheter is marked at 5, 10, 15 and 20cm from the distal end to fnd out the depth of placement.
- Proximal end of the catheter is ftted with Taper shape connector for easy connection to drainage Bottle/bag.

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 10; Master Carton of 100.
- Gross weight: 8.500kg
- Master Carton size: 560 x 400 x 260mm

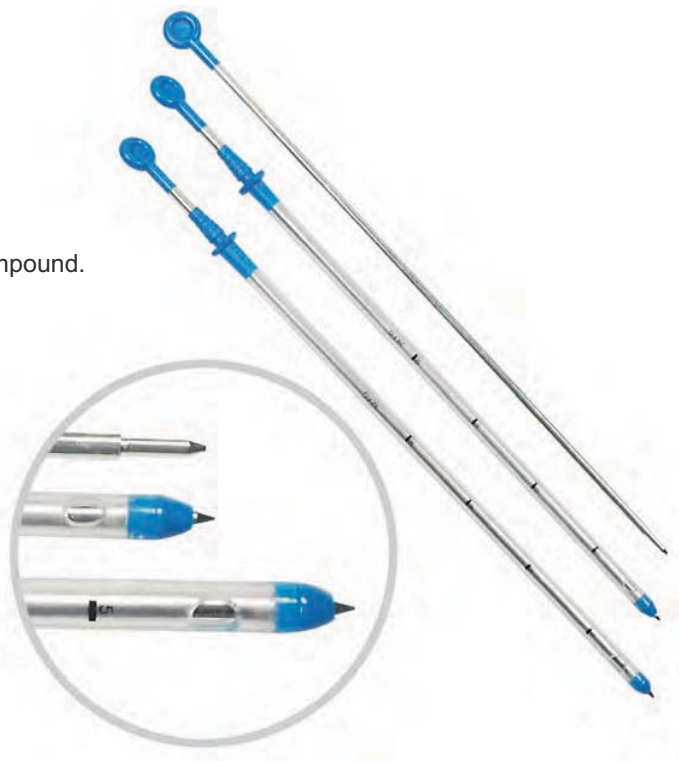

SIZES

Size	8	10	12	14	16	20	24	28	32	36
Length	230mm	230mm	230mm	230mm	230mm	390mm	390mm	390mm	390mm	390mm

THORACIC DRAINAGE CATHETER CHEST DRAINAGE CATHETER | MW 1002

- Thoracic Drainage Catheter is used for post operative drainage after cardio thoracic surgery.
- Manufactured from non toxic ,medical grade PVC compound.
- Smooth and round open distal end.
- Catheter is provided with large smooth eyes for maximum drainage.
- Catheter is marked at every 2 cm from the last eye to ascertain the depth of placement.
- Proximal end of the catheter is ftted with Taper tongue connector to providebetter forceps grip and smooth penetration
- Also available in right angle.

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 20; Master Carton of 200
- Gross weight: 11.500kg
- Master Carton size: 680 x 305 x 650mm



SIZES

Size	12	14	16	20	24	28	32	36
Length	450mm	450mm	450mm	450mm	450mm	450mm	450mm	450mm

STERI-ADK

ABDOMINAL DRAINAGE KIT | MW 1011

- Soft and smooth catheter with large atraumatic eyes for efficient drainage
- Radio-opaque line provided on the catheter for x-ray visualization
- Smooth and round distal end
- Catheter is marked at every 2cm from the cast eye to ascertain the depth of placement
- Proximal end of the catheter is provided with universal taper connector for easy connection to drainage bag.
- 2000ml bag manufactured from extra-strong sheeting
- Specially designed handle hold the tube upright
- Bag is graduated with easy to read scale

Packing

- Box of 5; Master Carton of 50pcs.
- Gross weight: 9.700kg
- Master Carton size: 570 x 320 x 340mm



STERI-DRAIN

CHEST BAG/WATER SEALED DRAINAGE SYSTEM | MW 1012

- Manufactured from non-toxic, non-irritant PVC
- Under Water sealed drainage system for collection of wound and empyema
- PVC Drainage bag with 1000ml capacity
- Disposable, Sterile ready for use.

Packing

- Individually packed in poly bag
- Master Carton size: 570 x 320 x 340mm
- Box of 5; Master Carton of 50pcs.
- Gross weight: 8.400kg



STERI-VAC SET

CLOSED WOUND SUCTION UNIT | MW 1003

- Steri-vac set is an effective device for closed wound drainage under negative pressure post operatively.
- Radio opaque drainage catheters are having multi perforation for efficient drainage.
- Suitable to offer surgeons/doctors to use one or two catheters simultaneously.
- Flexible Bottle chamber is easy to operate by one person to activate the suction.
- Available in different sizes with matching Needle.

Available With -

- Bottle Chamber with Hanger.
- Connecting tube with Matching size Connectors.
- Needle with perforated Catheter.
- One extra separate perforated Catheter.

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 1; Master Carton of 20.
- Gross weight: 7.200kg
- Master Carton size: 640 x 275 x 460mm



SIZES					
FG	10	12	14	16	18
Colour Code	Black	White	Green	Orange	Red

Chamber Capacity 800 ml and 500 ml

SURGICAL TAPES/SPORTS TAPES | MW 100

- Made from 100% cotton fabric for dressing and general medical applications.
- Coated with custom made non allergenic, non irritating adhesive.
- Adhesive is soft on the skin and does not pull out hair on removal.
- Thermo sensitive adhesive activate better when it comes in contact with warm human skin

Packing

- Packed in metal caps with plastic bobbin or paper core and card board in the form of rolls.
- Gross weight: 38.150kg
- Master Carton size: 432 x 482 x 648mm



SIZES & PACKING			
Width	Length	Pkg. (Per Box)	Pkg. (Master Carton)
12.5mm	1mtr.	45pcs.	36 Boxes
2.5mm	1mtr.	30pcs.	32 Boxes
2.5mm	5mtr.	30pcs.	6 Boxes
5.0mm	5mtr.	30pcs.	10 Boxes
7.5mm	5mtr.	10pcs.	10 Boxes
10.0mm	5mtr.	10pcs.	6 Boxes
30.0mm	10mtr.	25pcs.	12 Packs



SURGICAL PAD/GAUZE & COTTON TISSUE | MW 100

- Fine absorbent cotton covered with a layer of low limit gauze fabric.
- Meets BP/USP standards.

Packing

- Sterile, individually packed in peelable pouch pack.
- Carton of 500pcs.
- Gross weight: 7.320kg
- Master Carton size: 280 x 482 x 788mm

GAUZE SWABS | MW 102

- Made from fine surgical gauze/non woven fabric.
- Meets IP/BP standards.
- Optical brightener free.
- Free from foreign matter

Packing

- Available in Sterile/Non Sterile pack.
- Individually packed or available in combi pack of 2 or 4pc or 100pcs.
- Cartoon of 500pcs.
- Gross weight: 4.510kg
- Master Carton size: 584 x 381 x 254mm

SIZES
5cm x 5cm, 7.5 cm x 7.5 cm & 10cm x 10cm.
Option of 4 ply, 8ply, 12ply or 16ply.

SIZES
6cm X 6cm, 10cm X 10cm & 10cm X 20cm



STERIPORE | MW 106

- Steripore Surgical Tape is backed by a Thin highly permeable micro porous non woven fabric.
- Uniformly coated with hypo allergenic adhesive.
- Safe to use on sensitive part of the body.
- Suitable for infants and people with allergies.
- Peel off easily without pulling out hair or causing pain.
- High resistance to water and vapour and therefore will remain in position when immersed in water.

Packing.

- Carton of 40 Boxes
- Gross weight: 13.100kg
- Master Carton size: 535 x 320 x 270mm



SIZES & PACKING

Width	Length	Pkg. (Per Box)	Pkg. (Master Carton)
12.5mm	6mtr/9mtr.	24pcs.	40 Boxes
25.0mm	6mtr/9mtr.	12pcs.	40 Boxes
50.0mm	6mtr/9mtr.	6pcs.	40 Boxes
75.0mm	6mtr/9mtr.	4pcs.	40 Boxes



SIZES

10, 11, 12, 15, 20, 21, 22, 23, 24

SURGICAL BLADE WITH HANDLE MW 1006

- Disposable
- Ready to use, Sterile
- Fitted on handle No. 3 or 4

Packing

- Pack of 20. Carton of 200pcs.
- Gross weight: 1.780kg
- Master Carton size: 152 x 229 x 254mm

ABLECAST

ORTHOPAEDIC CASTING BANDAGE TAPE | MW 106

- ABLECAST has no harmful ingredients that could result in itching and a skin irritation unlike fiberglass tapes.
- Lightweight and well-ventilated ABLECAST can prevent terrible odors and infections.
- A good x-ray transparency makes it easier to take correct diagnosis

Packing.

- Box of 10. Carton of 60pcs.
- Gross weight: 16.300kg
- Master Carton size: 540 x 430 x 310mm

SIZES

Length 5cm, 7.5cm, 10.1cm, 12.5cm in 3.6mts



**SURGEON CAP AND NURSES CAP | MW 1204**

- Surgeon Caps/Nurses Cap are used to protect contamination from the head/hair.
- Surgeon Caps/ Nurses Cap are made from Non woven fabric.
- Provided with stitched interlocking.
- Elastic strip stitched on the back of the cap for proper grip on the forehead.
- Available in Blue /Green colour.

Packing.

- Individually packed .
- Box of 100; Carton of 1800pcs.
- Gross weight: 9.510kg
- Master Carton size: 660 x 432 x 508mm

SHOE COVER | MW 1210

- Shoe Cover is used as a covering of foot which serves as a protective barrier.
- Made on fully automatic ultrasonic sealing machine.
- Shoe Cover is made from non-woven fabric/polyethylene
- Elastic strip for proper grip.

Packing.

- Individually packed in polybag
- Carton of 2000 pairs. Pack of 25 pairs.
- Gross weight: 15.4kg
- Master Carton size: 737 x 584 x 432mm

SURGEON FACE MASK | MW 1208

- Surgeons Face Mask is used as a covering for the face that serves as a protective barrier.
- Surgeons Mask is made from Non woven fabric.
- Available in two layers and three layers configuration.
- High bacterial filtrations efficiency Ultrasonically welded.
- Available in Tie bands and Elastic Strap for proper fixing.
- Provided with a nose clip to prevent fogging of spectacles and to maximize filtration efficiency.

Packing.


- Individually packed
- Box of 100; Carton of 4000.
- Gross weight: 19.570kg
- Master Carton size: 660 x 432 x 508mm

SURGICAL GOWN | MW 1202

- Surgical gown are made from Non woven fabric to provide softness and comfort during use.
- Provided with stitched interlocking.
- Soft elastic cuff
- Tie straps are provided for better fitting and complete body wrap
- Light in weight

Packing.

- Individually packed
- Box of 50.
- Master Box of 100kg.
- Gross weight: 15.00kg
- Master Carton size: 533 x 533 x 558mm

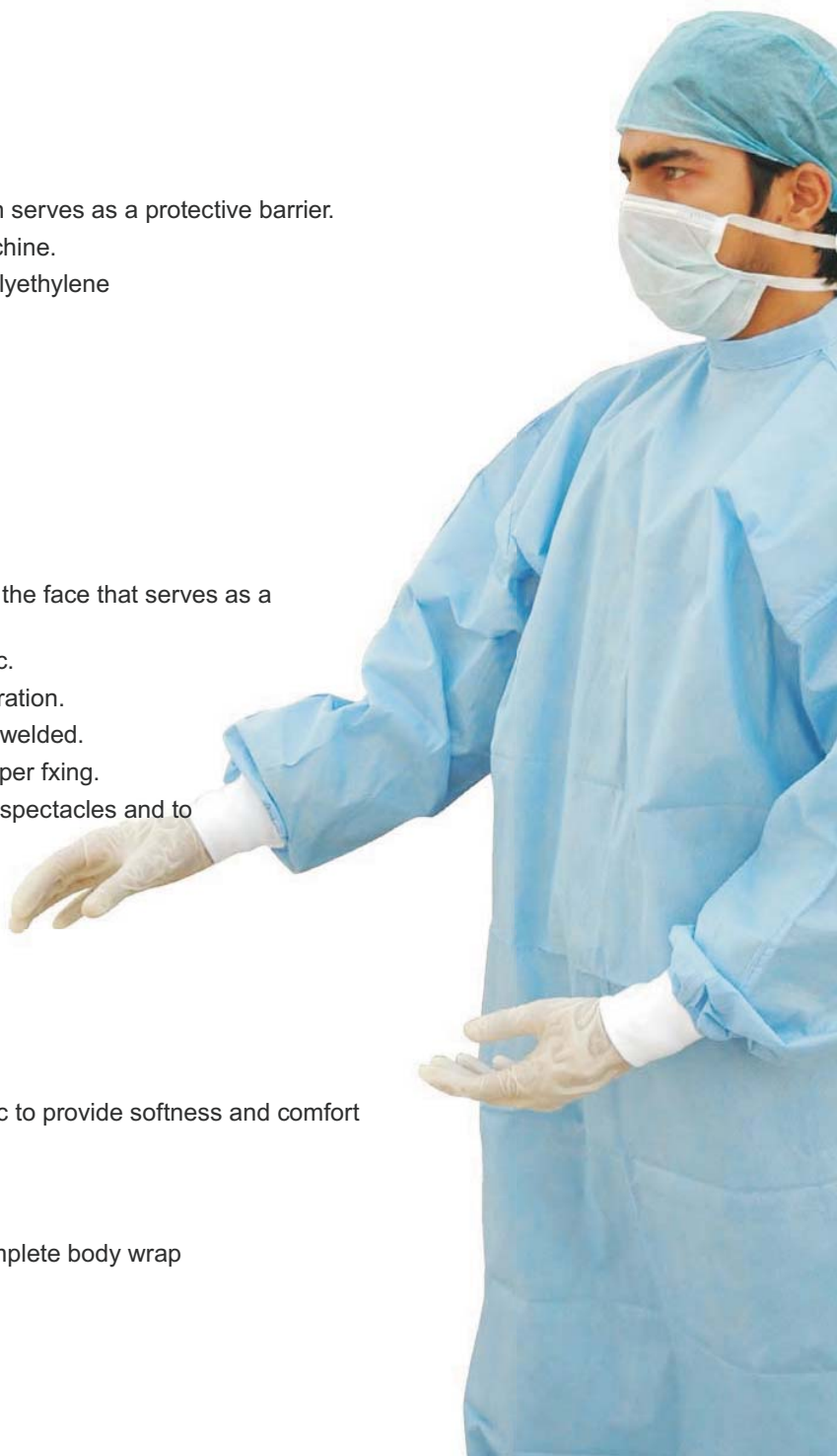



H.I.V KIT | MW 1013

- Disposable
- Kit consist of

- Proto Gown	-1Pc.
- Hood Cap	-1Pc.
- Surgeon's Cap (with Bacterial Filter)	-1Pc.
- Shoe Legging (Knee Length)	-1Pair
- Latex Gloves	-1Pair.
- PE Sheet	-1Pc.
- Goggles	-1Pc.
- Carry Bag	-1Pc.
- Face Mask	-1Pc.
- We can also make tailor made sterile disposable kits as per the Customer's requirement such as
 - Gynecology & Obesteries
 - General Surgery
 - Orthopaedic
 - Urology Kit

Packing

- Sterile, individually packed in peelable pouch pack.
- Box of 10.

PATHOLOGY LAB PRODUCTS

URINE/STOOL SPECIMEN CONTAINER

- Urine/Stool Specimen Container is used for the Collection of Urine /Stool Sample for the laboratory examination.
- Made from medical grade Translucent Polypropylene or clear polystyrene

Packing

- Sterile, individually packed in Tear pouch pack.

URINE SPECIMEN CONTAINER | MW 901

- Urine Specimen Container provided with easy to read graduation helps to determine the Volume.
- Provided with airtight Lid for safe handling of sample
- Available in 30ml and 70ml capacity.
- Pack of 25. Carton of 600.
- Gross weight: 9.230kg
- Master Carton size: 660 x 432 x 508mm



STOOL SPECIMEN CONTAINER | MW 902

- Provided with airtight spoon type Lid for safe lifting & handling of sample.
- Pack of 25.



BLOOD LANCET | MW 900

- Blood Lancet used as a manual blade to make a limited skin incision as for collection of Blood.
- Made from medical grade Stainless Steel and polymer.
- Ultrasonic cleaned.
- Super sharp tip for the instant incision.

- Master Carton size: 432 x 432 x 432mm

Packaging

- Sterile, individually packed in peelable pouch pack.
- Pack of 200. Carton of 20000pcs.
- Gross weight: 13.00kg



... enriching life

V. S. BIOSCIENCES LLP

A-204, Neelam Centre,
Hind Cycle Road, Worli,
Mumbai - 400 030 (INDIA)

Tel.: (+91 22) 2498 8777
Tel.: (+91 22) 4333 3900
Fax: (+91 22) 6660 8779